

PORTER

# Silhouette

LOW PROFILE NASAL MASK

NEW &  
IMPROVED



**Enhanced Features  
New Lower Price**

**PORTER**

*The Trusted Name In Nitrous*

[PorterInstrument.com/Silhouette](http://PorterInstrument.com/Silhouette)



Robert Guyette,  
D.M.D., M.D.

## Designed By An Oral Surgeon To Make Nitrous Oxide Easier To Use

The Porter SILHOUETTE nasal mask is the brainchild of Robert Guyette, D.M.D., M.D., an Oral and Maxillofacial Surgeon in Scottsdale, AZ. Dr. Guyette recognized that the available nitrous oxide masks made it difficult to access the oral cavity during procedures which resulted in him avoiding using nitrous oxide for his patients.

His solution was to design an anatomic, low profile mask that fits comfortably over the patient's nose while also allowing his patients full range of head motion – and SILHOUETTE was born!

Now, after years of successful application in dental offices around the world, SILHOUETTE has been redesigned with a number of enhancements, making it even better – along with a new lower price.

## Description

Like its predecessor, Porter's new and improved SILHOUETTE nasal mask is a translucent, single-use nasal mask and breathing circuit designed with lightweight tubing and a comfortable, conforming mask. It provides unobstructed access to the oral cavity and full range of motion for the patient.



Traditional nasal mask




Porter SILHOUETTE nasal mask

## Features & Benefits

The new SILHOUETTE mask is narrower, flexed in, and combined with the new support bars, pulls the mask in against the nose. Other advantages over traditional nasal hoods include:

- Conforming shape and fit without adhesive
- Enhanced scavenging efficiency
- Colored masks for easy size identification
- Disposable mask and circuit for improved infection control
- Four size options
- New nasal barb shape and size for improved flow and comfort
- Mild mint scent
- Capnography port in each mask
- No tubing behind headrest
- Protective eyewear fits over mask

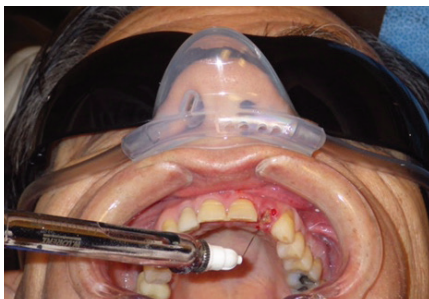

SILHOUETTE allows unobstructed access to the oral cavity and full range of motion.

# SILHOUETTE Low Profile Nasal Mask And Breathing Circuit

Fits securely without adhesive

New and improved mask design

Built-in capnography port



Single-use mask and tubing for improved infection control and efficient room turnover



## A Comfortable Fit For Patients



The new and improved SILHOUETTE mask is designed to fit comfortably and securely over the bridge of the nose without adhesive.



## Scavenging Efficiency

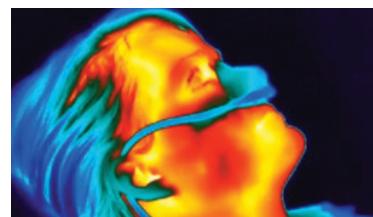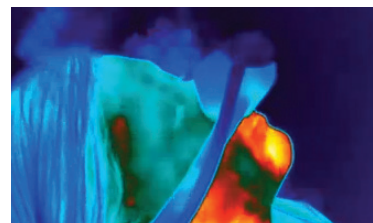

Visit our website to see videos of SILHOUETTE's scavenging efficiency.



Mask redesigned to pull in on the nose and fit securely. Cannula shape and size redesigned to improve flow as well as a smaller size for small and pediatric masks.



Patients can wear protective eyewear and maintain an excellent seal of the mask.



Tubing is in front of the patient (versus behind headrest) which allows for full range of head motion. Patients can turn their head, sit up, and rinse without the mask becoming dislodged.

**VISIT OUR WEBSITE FOR VIDEOS, PHOTOS AND FAQs – [WWW.PORTERINSTRUMENT.COM/SILHOUETTE](http://WWW.PORTERINSTRUMENT.COM/SILHOUETTE)**

## How To Get Started

It's quick and easy to start using Porter's SILHOUETTE nitrous oxide nasal masks. SILHOUETTE can easily be retrofitted to all dental flowmeters. Get started today by ordering a Retrofit Package and Variety Pack!



## PORTER SILHOUETTE PART NUMBERS

PART	DESCRIPTION
SIL2-RETRO-PKG	Flowmeter Retrofit Package (order one for each flowmeter) Includes: Cannula Adapter, Breathing Bag Cap, 5.5' Fresh Gas/ Vacuum Tubing, 4 Color-Coded Sizer Masks, and Clip
SIL2-VAR-4X3	Variety 12 Pack: 3 Silhouette Circuits of each size - Pediatric, Adult Small, Medium, Large
SIL2-PEDO-12 or SIL2-PEDO-24	Pediatric 12 Pack or Pediatric 24 Pack
SIL2-SM-12 or SIL2-SM-24	Small 12 Pack or Small 24 Pack
SIL2-MED-12 or SIL2-MED-24	Medium 12 Pack or Medium 24 Pack
SIL2-LG-12 or SIL2-LG-24	Large 12 Pack or Large 24 Pack
SIL-SIZER-4	SILHOUETTE Sizer Masks - Color Coded, 1 For Each Size
SIL2-CAPNO	SILHOUETTE Capnography Barbs (25)
SIL-VAC	SILHOUETTE Vacuum Adjuster



SILHOUETTE easily retrofits to any flowmeter with SIL2-RETRO-PKG. See retrofit videos on our website.

Order SIL2-CAPNO for capnography monitoring.

SILHOUETTE is only available through authorized Porter distributors.

### Contact Porter for more information:

P: 215-723-4000 Email: porternitrous@parker.com



[www.porterinstrument.com/silhouette](http://www.porterinstrument.com/silhouette)